





SIMPLE CONTROL

Chart & Projector ULC-900

The **NEW** ULC-900 visual acuity chart system, eye level and distance can be adjusted by simple control, and it gives more convenient environment for users.

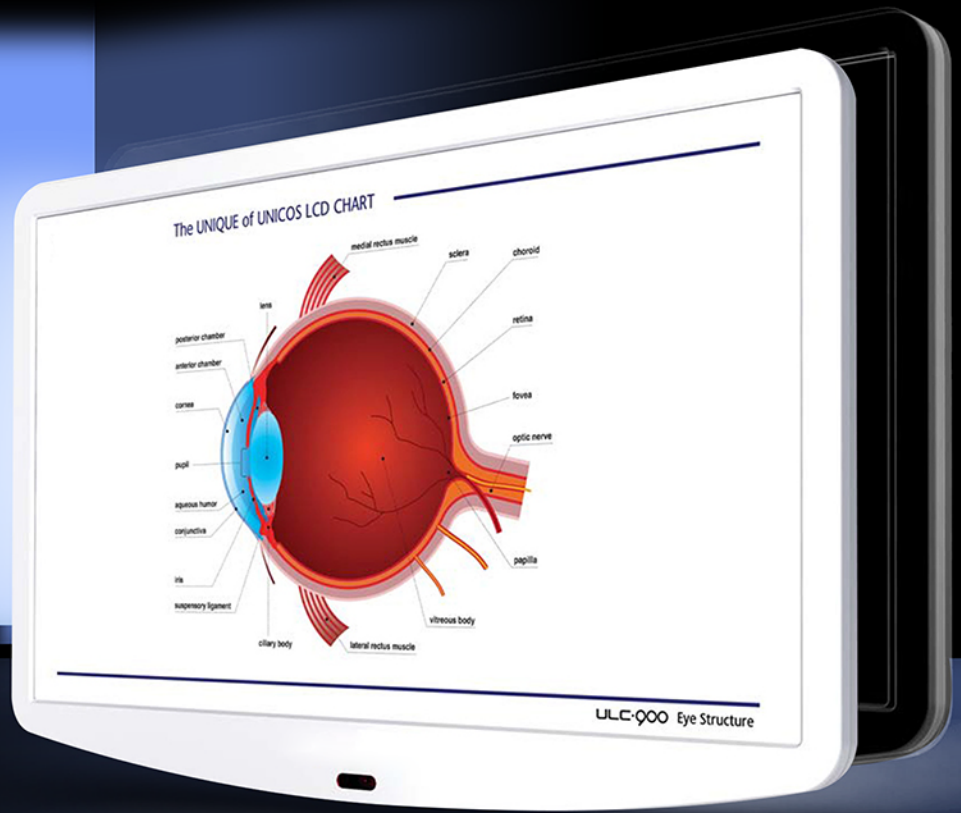

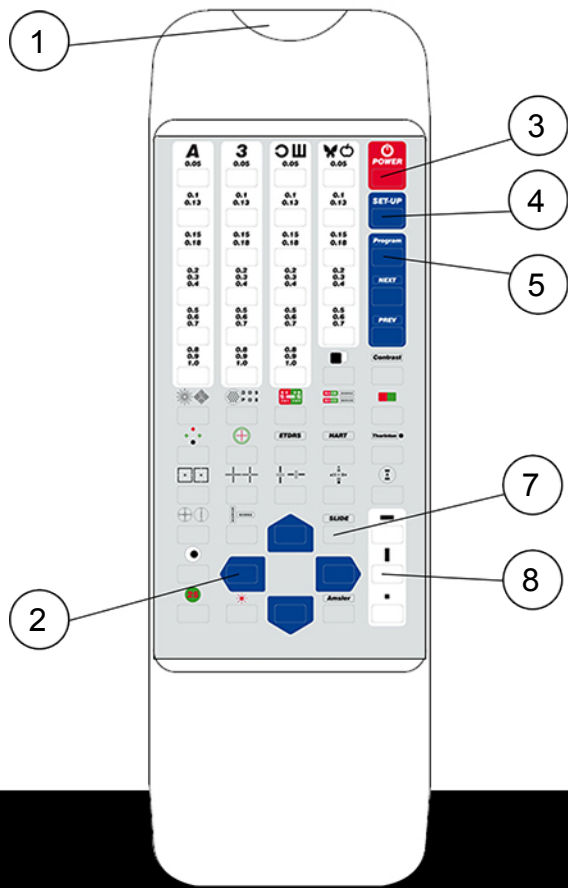


**Chart & Projector
ULC-900**



REMOTE CONTROL

ULC-900 features a remote control designed for prompt and convenient tests. Easy image addition at slide mode from USB flash drive.

- 1. Transmitter
- 2. Arrows
- 3. Power Button
- 4. Set-Up Button
- 5. Program Button
- 6. Slide Button
- 7. Mask Button
- 8. Serial Cable Port
- 9. Power Switch
- 10. USB Port
- 11. Adapter Connector
- 12. Wall Mount Bracket (Optional)

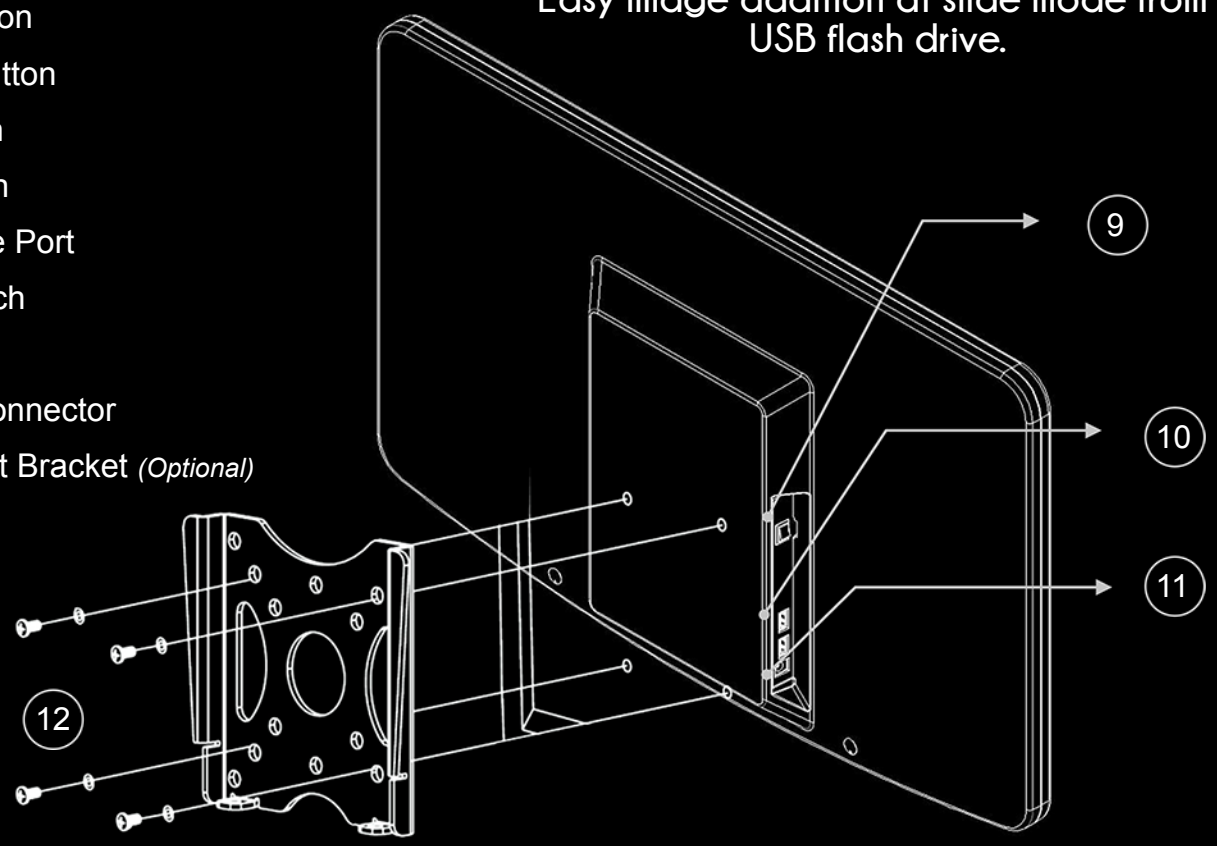

ULC-900 | Chart & Projector


24 Inch TFT Monitor

High-resolution display for wider and clearer view.

Distance: 5.00 (m)
Contrast: 5 (%)
Unit of Index: Decimal
 Feet
 Meter
ETDRS Chart: TYPE 1
 TYPE 2
 Landolt
RG / Polar: Polar
Cross Type: Type 1
 Type 2
Mask Position: Original
 Center
Mirror: ON
 OFF
RS-232 Protocol: B Type
Remote Channel: (Ch)
CH: (Ch)
LED
LCD Brightness
Power Save
Set Default
Reload
S/W Update
Update

High Performance Polar Filters

- Keep charts clear between the range of 1.5M ~ 7.5M test distance by high performance polar filters. Heterophoria Heterotropia / Stereopsis tests with polar separation charts.

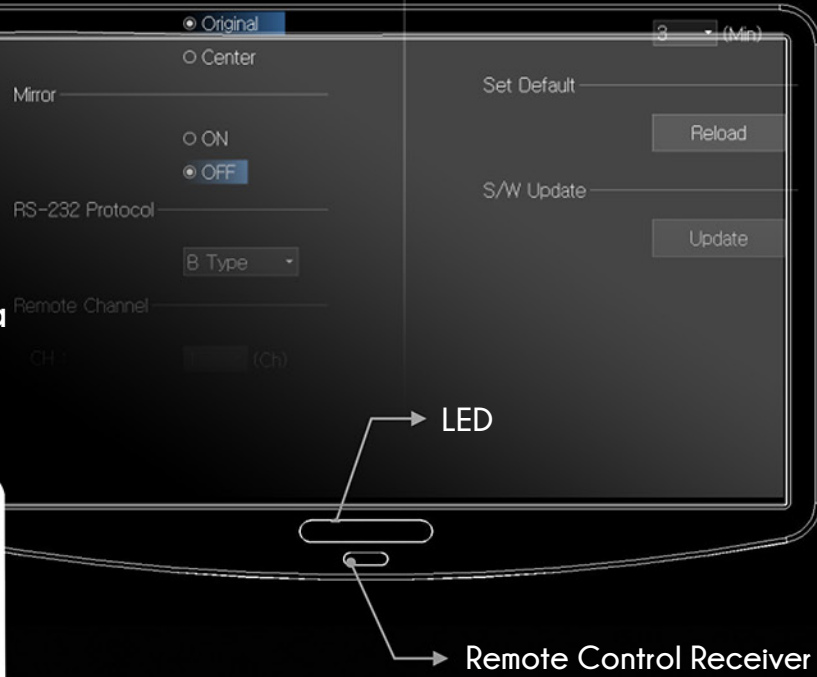

Installation Distance

May control installation space from 1.5M up to 7.5M by 0.25m. May use of Mirror functions at the small place.



GENERAL CHARTS

ULC-900

Chart & Projector


Back 0.05 ~ 2.00 charts by test distance up (Children 0.1~1.0) Horizontal, Vertical, Single character mask, Center indication Red / Green filters / Random charts.



Contrast



HART

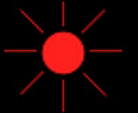

ETDRS



D 9
P 8
R 6



Thorinton



Configuration

- Chart
- Red / Green
- Program
- Image

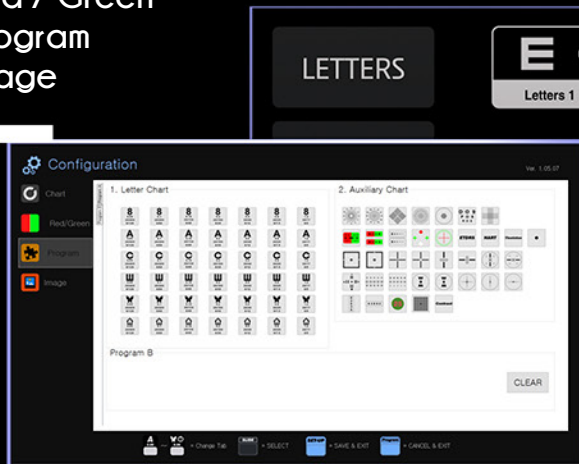


Chart Test Options

- Binocular balance test, Astigmatism & Stereoscopic photos

White-black contrast

- Making use of testing Lasik operated patient.

Filter level configuration


- Correct filter levels with MENU function may check corrected levels on monitor.

Numbers of different charts

- Landolt Ring, Snellen, Child chart, Cross chart and more.

Color-blindness test












- Preset seven type's charts / Random selection function.


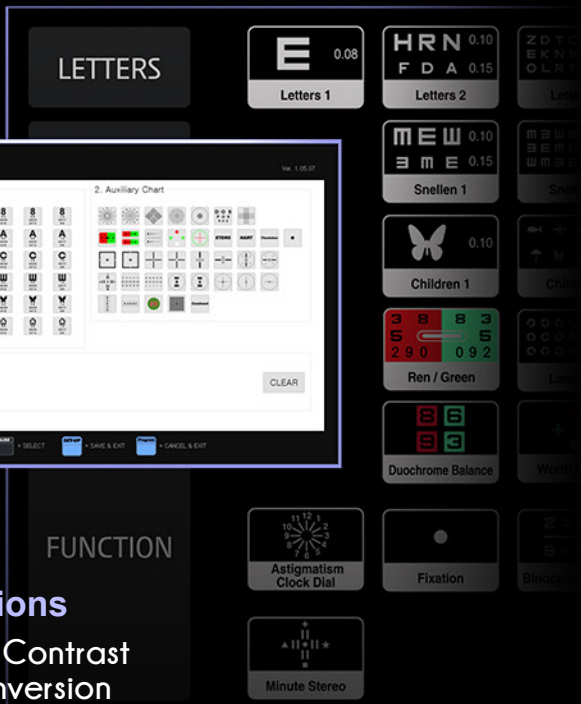


Special Functions

- Mirror mode / Contrast adjustment / Inversion



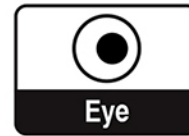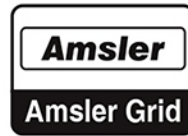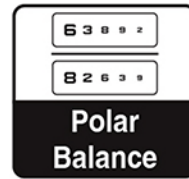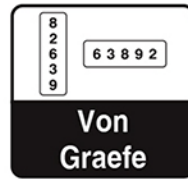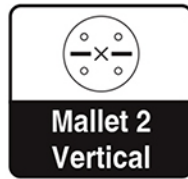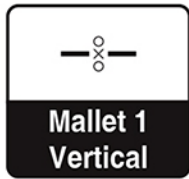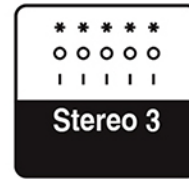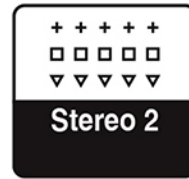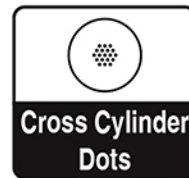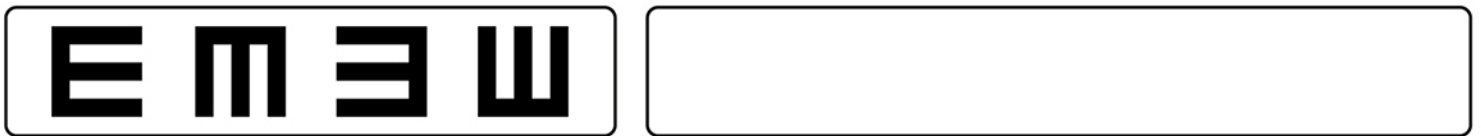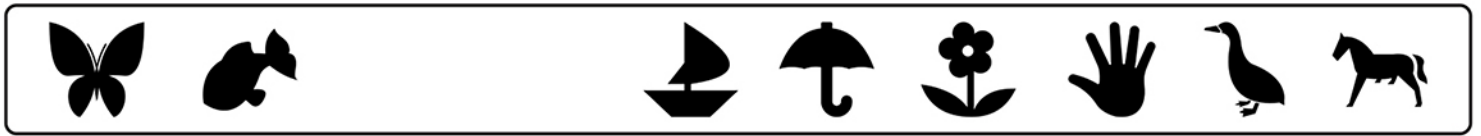

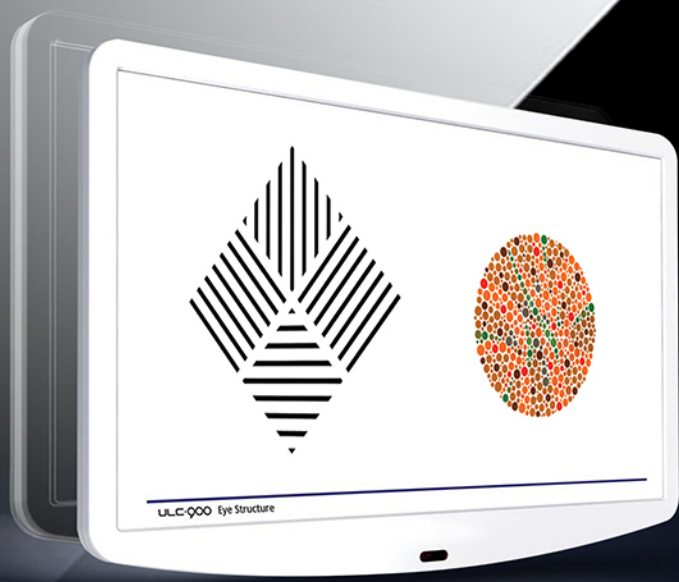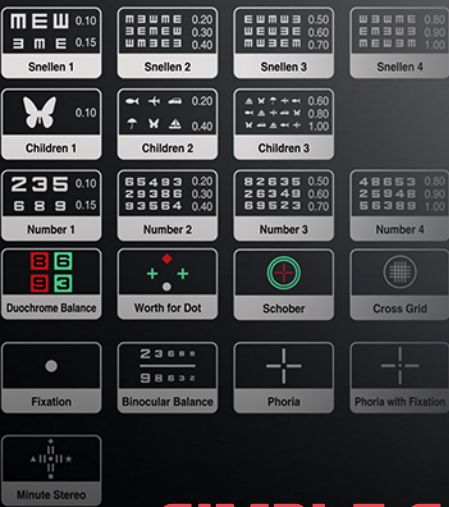

SIMPLE CONTROL

ULC-900

Chart & Projector

Specifications

LCD Type	24Inch TFT, LED Backlight
Resolution	1920 (H) X 1080 (V) Pixels
Display Area	531.36 (H) X 298.89 (V) mm
Dimension / Weight	575(W) X 370(H) X 43(D) mm, 4Kg
RC size / Weight	64(W) X 195(H) X 20(D) mm, 110g
Power Supply	DC 5V / 5A Adapter
Power Consumption	Max. 15W
Chart	Landolt Ring / Snellen Numbers / Letters / Korean Children
Masking	Single Character, Horizontal Line, Vertical Line Normal, Center
Filter	Red / Green Filter
Visual Functions	Red / Green, Polar separation, 3D polar, Strabismus, 3D Red / Green, Color blindness
Distance	1.5M ~ 7.0M (Changes into 0.25meter unit)

Accessories

- Cradle - Wall (Standard) / Floor (Optional) / Table (Optional)
- Special glasses - Polarizing glasses (Optional),
- Red/green glasses (Optional) /Remote Controller (Batteries), Adapter,
- Serial cable, Power cable, Cleaning cloth

Designs and details can be changed without prior notice for its improvements.

