

Inami MANUAL












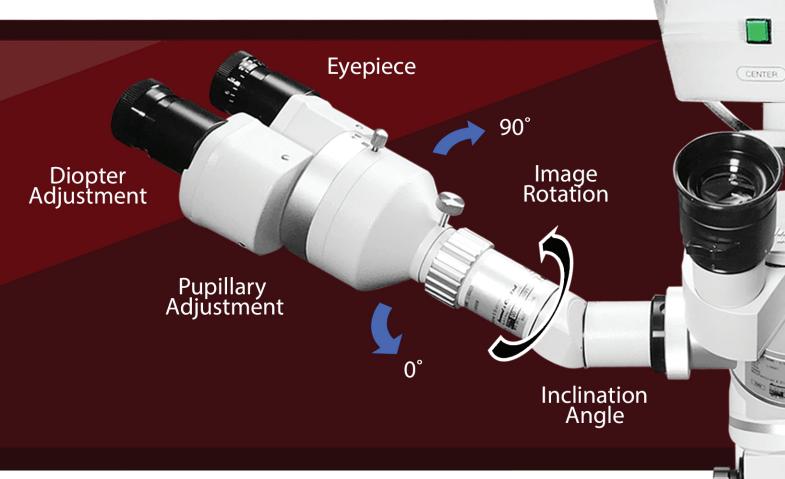



Specification

Eyepiece : 12.5 X

Diopter adjustment : $+5D\sim-6D$

Pupillary adjustment: 48~78 mm

Image Rotation : $\pm 100^{\circ}$ Inclination : $0 \sim 90^{\circ}$

Mirror Cleaning

Remove dust accumulation with a camel's hair brush.

After brushing, finger marks can be removed, if necessary, with isopropyl alcohol and a cotton swab.

Dry with a facial tissue.



Microscopes Accessories



This ergonomic design microscope features an optional C mount system, beam splitter, assistant binocular microscope, and video USB camera.





Included














Inami microscopes offers an optional USB digital camera, with a unique software where you can store and create your own patient data base to capture high quality digital photos and videos.

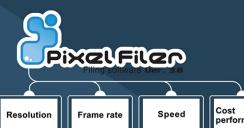







Assistant Binocular Microscope

The optional high quality assistant binocular coaxial microscope offers an apochromatic optics which allows the observer to magnify and visualize the red flex with easy and control.



Beam Splitter Devices

All operation microscopes features an optional beam splitter that can be easily assembled and it allows you to attached many other Inami optional accessories.



Video Camera Adaptor C-Mount F=70mm

This optional video digital camera adapter is designed for a easy and quick assembly access while adding the USB digital camera.